The logo for Wellbox X2 features the word "wellbox" in a dark grey, lowercase sans-serif font. To its right is a stylized "X" composed of two overlapping chevron shapes: a red one pointing right and a dark grey one pointing left. A red "2" is positioned to the upper right of the "X".

wellbox X²

The "MADE IN FRANCE" logo consists of the words "MADE", "IN", and "FRANCE" stacked vertically in a white, uppercase sans-serif font. To the left of the text is a small graphic of the French flag's vertical stripes: blue, white, and red.

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THANK YOU
FOR CHOOSING
THIS DEVICE
WELLBOX® 2

We're delighted to have you as one of our customers!

The LPG® Group, world leader in cellular stimulation with its professional technologies, has developed the Wellbox® 2 device so that you can benefit from its expertise at home.

Wellbox® 2 uses different massage heads to care for your body and face.

It can be used to:

- Slim the figure
- Smooth cellulite
- Firm the skin
- Regain light legs
- Reshape the contours of the face
- Reduce wrinkles
- Restore radiance to the complexion

Compact and intuitive, the Wellbox® 2 device includes 80 care and programmes via its touch screen to meet every beauty and well-being objective.

The Wellbox® 2 device is designed for home use by the general public; it is not intended for professional use.

USER'S GUIDE

To reduce environmental impact while maintaining safety and ease of use, the full user guide is available in electronic format at www.wellbox.com


Be sure to read this guide before using your device.

If you have any doubts about the operation or use of your device, or if you would like to receive these instructions in hard-copy format, please contact Wellbox® Customer Service:

+33(0)4 75 78 69 00 (toll-free in Europe only)

PRODUCT WARRANTY

ACTIVATE YOUR 2-YEAR GUARANTEE ON THE WEBSITE: WWW.WELLBOX.COM OR BY SCANNING THE QR CODE OPPOSITE USING YOUR PHONE'S CAMERA.

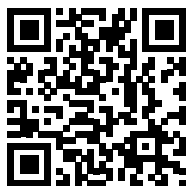PACKAGE CONTENTS

The Wellbox®2 device comes with all the equipment you need for body and facial care.

- 3 ROLL care heads (1,2,3)
- 2 LIFT care heads (4,5)
- A connection hose for connecting the care heads
- A power supply unit
- Two disassembly tools for the LIFT 4 and 5 heads
- Replacement filters
- Replacement flaps for the LIFT heads
- A quick user guide

IMPORTANT: Keep the original Wellbox® 2 packaging for the duration of the warranty. Any returns during this period must include the original packaging.

CONTRAINDICATIONS

Do not use the Wellbox® 2 device if you have any medical contraindications, or if you have not read this chapter carefully and reviewed the precautions in this guide. Consult your doctor if you have any doubts about the compatibility of your health with the use of the Wellbox® 2 device.

The Wellbox® 2 device can be used on all healthy external parts of the body, with the exception of the eyes, ears, glottis and genitals.

The Wellbox® 2 device should not be used on mucous membranes, nipples, the spine or on any suspicious moles/beauty spots without medical advice.

If you are pregnant, do not use the device on your stomach/abdomen or lower back. Some medical conditions may require the advice of your doctor.

Do not use the Wellbox® 2 device if you have any of the following:

- Epilepsy problems
- Blood clotting disorders
- Heart problems
- Pain or redness
- Tingling in the limbs
- Phlebitis, varicose veins or if you have undergone an operation on your veins in the previous two months.
- Fever, infection, inflammation, acute dermatitis, open wounds, recent scars or any serious and/or progressive illness.
- In case of any prolonged use of medications (anti-inflammatories, anticoagulants, antibiotics and antihistamines), consult your doctor before using this device.

Discontinue care if you experience any unexplained pain in your calves.

If you experience any side effects related to the Wellbox® 2 care during or after use, consult your doctor and contact Wellbox® Customer Service.

If you have any doubts about the compatibility of your health with the use of the Wellbox® 2 device, consult your doctor.

PRECAUTIONS

- Do not use the Wellbox® 2 device on children under the age of 7. Do not allow a child to use the Wellbox®2 device alone. Remember to reduce the intensity when treating a child.
- Close supervision is necessary when the device is used by, on or near children, invalids or disabled people.
- This device may be used by children aged 11 and over and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instructions to use this device safely, and provided they understand the hazards involved.
- Do not allow children to play with the device.
- Cleaning and maintenance must not be carried out by unsupervised children.
- **Using the Wellbox®2 device should never be painful or leave marks; any redness should not last more than a few seconds or leave traces.**
- Never use the Wellbox®2 device without its protective filter.
- Use this device only as described in this guide. Use only the care heads and accessories supplied by LPG SYSTEMS.
- Do not use detergents, solvents, chemicals or abrasive sponges to clean the Wellbox® 2 device. Use wipes soaked with a bactericide and fungicide solution.
- Never use the Wellbox® 2 device on jewellery or piercings.
- Do not operate the device at altitudes above 2,000 metres.
- Under normal conditions of use, the device should be stopped for 30 minutes between two 30-minute sessions.
- Do not use the device in an oxygen-rich environment.
- The device should always be positioned so that the mains plug is easily accessible.
- **The device must always be used on clean, healthy skin (no creams, oils, etc. before use).**

Using the Wellbox® 2 device regularly will help you feel better; we recommend not interrupting its use for too long. Failing to follow the instructions and recommendations may prevent you from obtaining the desired results. Do not use the Wellbox® 2 device for longer than the recommended time. **Results may vary depending on your age, your lifestyle (sporting activity, diet, etc.) and your morphology.**

IMPORTANT SAFETY INFORMATION

.....

Before installation and use, you are strongly advised to read the instructions in this chapter.

DANGER - TO REDUCE THE RISK OF ELECTRIC SHOCK:

- Check that the mains voltage is consistent with the range indicated for the device's power supply (100-240 VAC). The device must be connected to a wall socket that complies with your country's specific standards. Avoid connecting the Wellbox® 2 device to the mains via a power strip.
Disconnect the Wellbox® 2 device from the power supply after each use or before cleaning.
- Do not use the Wellbox® 2 device while taking a bath or shower. Do not place the Wellbox® 2 device on a damp or unstable surface or near to a bathtub, shower or sink. Use the Wellbox® 2 device in a dry place.
Avoid splashing any water on the Wellbox® 2 device. Do not touch any device that has fallen into water; disconnect the power supply immediately.
Do not immerse the hose in water. Do not use the Wellbox® 2 device outdoors.
- Do not drop or insert objects into the openings of the device.
- Do not allow any solid object or liquid to be sucked into the Wellbox® 2 device.
- The device must always be used only with the power supply unit provided.
- The device must be powered only with the safety extra-low voltage corresponding to the marking on the device.

IMPORTANT SAFETY INFORMATION

(CONTINUED)

.....

WARNING - TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK OR INJURY:

- Do not use the device under a blanket or similar item.
Excessive heat could occur and cause fire, electric shock or injury to persons.
- Never use the device if it has a damaged cable or plug, if it is not working properly, if it has been dropped or damaged, or if it has fallen into water. In these cases, contact Wellbox® Customer Service.
- Do not carry or move the device by its top cover, cables or hose, or use them as carrying handles.
Do not drop the Wellbox® 2 device.
- Do not shake or knock the device.
- Do not use the Wellbox® 2 device without first connecting the hose to the heads. Keep cables and hoses away from sources of heat.
Do not use the device when the ambient temperature is higher than 40°C.
- Do not use the device after using aerosol products (sprays) or when oxygen is being administered.
- Do not leave the device unattended while plugged in.
To switch off the device, press the on/off button, then disconnect the AC/DC power supply unit from the mains socket.
- Do not open or disassemble the device.
- The Wellbox® 2 device must be powered only with 24V DC safety extra-low voltage.

FIRST INSTALLATION

- 1 Open the device by pressing the push-button on the front.
- 2 Store the 5 care heads in their dedicated slots, recognisable by their markings.
- 3 Insert the hose with the raised arrow pointing upwards as shown in the photo.
- 4 Choose one of the care heads (LIFT or ROLL) and connect it by aligning the red arrows.
- 5 The Wellbox® 2 device is supplied with a power supply unit. Input 100-240 VAC, 50/60Hz - Output 24 VDC. Connect the power cable to the back of the Wellbox® 2 device. Then plug the power supply unit into a wall socket.
- 6 Once you have plugged in your device, the screen will turn on automatically. Now you can enjoy your device!

1 2 3

4 5 6

For optimum performance, we recommend that you leave your device plugged in for at least 2 hours:

- After the first use
- After any long absence of use (6 months)

DESCRIPTION OF THE TOUCH SCREEN

The Wellbox® 2 device is equipped with a touch screen with access to 5 menus:

MY ACTIVITY: find your current programmes, your care history, and track your progress.

TARGETED CARE: work on the area (body or face) of your choice depending on your objective (cellulite, slimming, radiance, firmness, etc.)

PROGRAMMES: follow a programme for several weeks (face, body or body & face)

SETTINGS: define your settings and get information on your device

M PROFILE: modify your profile or create a new one (maximum two profiles)

For step-by-step guidance, see the full user guide available on our website.

Watch video advice from the Wellbox® 2 expert by scanning the QR code opposite using the camera of your phone.

GUIDANCE FOR USE

BODY CARE

The **red lines** indicate the position of the rollers as they are moved.
The **arrows** indicate the direction of movement of the care head.


Two directions of movement are available:

ONE-WAY MOVEMENT:

The care head moves in the direction indicated and then returns to the starting point.
To be repeated throughout the care session.

THERE & BACK MOVEMENT:

When using the ROLL 1 head, press the centre of the trigger to reverse the direction of the rollers.
The care head must remain in contact with the skin throughout the care.

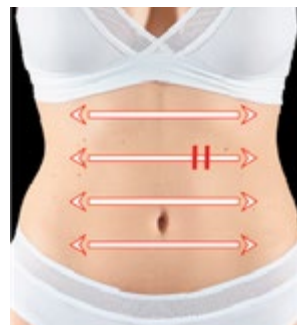Trigger

ROLL 1 HEAD


The number of movements will depend on the size of the area to be treated.

FACIAL CARE

The **red lines** indicate the position of the flaps during the movement.
They must always be perpendicular to the wrinkle.

For improved stimulation, **move your LIFT head in a stop-and-start movement, stopping for 4 seconds** at each point.

The number of movements will depend on the size of the area to be treated.

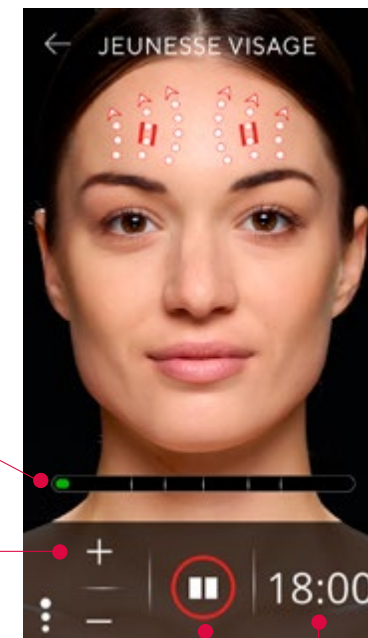The progress gauge represents the time elapsed for each stage.

Everything is pre-set, simply adjust the suction power by pressing the **"+" and "-" buttons** and follow the instructions on the display.

Press **PLAY** to start the care. You can pause it at any time by pressing the **PAUSE** button.

Total duration of the session

- The suction power should be adjusted according to the quality of your skin. Be careful not to pull on your skin. The care must **ON NO ACCOUNT** cause pinching or pain. Intensity is not synonymous with results.
- Regular use of your Wellbox® 2 device is the key to good results. We recommend starting with an “attack” phase, every day for anti-ageing/firmness, and every other day for slimming/anti-cellulite.
- For optimum effectiveness, it is important to have a healthy lifestyle, i.e. a balanced diet and regular physical activity (such as walking).
- Applying LPG® cosmetics after each session helps to optimize results. Once stimulated, the skin is more receptive to the active ingredients in the cosmetics.

IMPORTANT do not apply any cosmetics before the session, as this may damage your device.

Results can be more or less rapid depending on each user's condition (pre-menopause, hormonal influences, medical problems) and metabolism.

STORAGE

To close the Wellbox® 2 device after use:

- 1 Connect the ROLL 1 head to the hose.
- 2 Wind the hose into the purpose-built groove (2 turns).
- 3 Then place the ROLL 1 head in the space reserved for it next to the touch screen.

In case of any difficulty, we recommend removing the ROLL 1 head from the hose or removing the hose completely from the Wellbox® 2 device.

1

2

3

MAINTENANCE

In the interests of hygiene and to ensure proper operation of your device, please:

- 1 Clean the care heads after each use using wipes soaked with a bactericide and fungicide solution.

- 2 Replace the protection filter when indicated on the screen.

Watch the video advice on how to maintain your Wellbox® 2 by scanning the QR code opposite using your phone's camera.

WARNING

It is strictly prohibited to modify your device or perform any technical intervention on it without the authorisation of LPG SYSTEMS.

TECHNICAL SPECIFICATIONS

Dimensions L x W x H: 350 x 300 x 250 mm

Net weight: 7.5 kg

Power consumption: 70 Watts

Class 3 device, maximum duration of use per session: 30 minutes

Time required for the Wellbox® 2 device to warm up if stored at -25°C: 9 hours

Time required for the Wellbox® 2 device to cool down if stored at +70°C: 18 hours

Maximum surface temperature of the Roll 1 head: 42.5°C

Maximum surface temperature of the rollers on the Roll 1 head: 44.5°C

Protection rating _____ IP21

Electrical characteristics _____ 24V; 3.75A

Maximum regulated vacuum _____ 350mbar

Maximum altitude _____ 2.000 m

TECHNICAL CHARACTERISTICS OF THE POWER SUPPLY UNIT

Manufacturer _____ MeanWell

Ref. _____ GSM90B24-P1M

Electrical protection class _____ II

Protection rating _____ IP22

Electrical specifications _____ 100-240VAC 50/60Hz 1.5A

ENVIRONMENTAL CONDITIONS OF USE

Temperature _____ +5°C to +30°C

Relative humidity _____ 15% to 90%
(without condensation)

Atmospheric pressure _____ 700 hPa to 1.060 hPa
(for use in a normally ventilated room).

ENVIRONMENTAL CHARACTERISTICS OF TRANSPORT AND STORAGE

Temperature _____ -20°C to +70°C

Relative humidity _____ 10% to 90%
(without condensation)

Atmospheric pressure _____ 500 hPa to 1.060 hPa

IDENTIFICATION LABEL

"Type" identification label subject to change to take account of upgrades to the device and to the Wellbox® device type.
The label on the bottom of your Wellbox® device is the approved version.

WARNING
The manufacturer reserves the right to change the product specifications of the device.

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SIÈGE SOCIAL LPG SYSTEMS SAS
LPG SYSTEMS SAS HEAD OFFICE