# wellbox<sup>2</sup>

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USER GUIDE

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MANUALE D'USO

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# THANK YOU FOR CHOOSING THIS DEVICE **WELLBOX® 2**

.....

We're delighted to have you as one of our customers!

The LPG® Group, world leader in cellular stimulation with its professional technologies, has developed the Wellbox® 2 device so that you can benefit from its expertise at home.

Wellbox® 2 uses different massage heads to care for your body and face.

It can be used to:

- Slim the figure
- Smooth cellulite
- Firm the skin
- Regain light legs
- Reshape the contours of the face
- Reduce wrinkles
- Restore radiance to the complexion

Compact and intuitive, the **Wellbox® 2 device includes 80 treatments and programmes** via its touch screen to meet every beauty and well-being objective.

Be sure to read this guide before using your device.

If you have any doubts about the operation or use of your device, please contact Wellbox® Customer Service:  
00 800 699 355 269, Monday to Friday, 9am to 5pm.

**The Wellbox® 2 device is designed for home use by the general public; it is not intended for professional use.**

# CONTENTS

Wellbox® 2 warranty certificate ..... 31

Contents of packaging ..... 32

Contraindications ..... 33

Precautions for use ..... 34

Important safety instructions ..... 35

Description of the device ..... 37

Description of the touch screen ..... 38

Description of the care heads ..... 39

Initial installation ..... 41

Getting started ..... 42

Storage ..... 44

Maintenance ..... 46

Roll head maintenance ..... 47

Lift head maintenance ..... 48

Problems? What to do? ..... 49

Technical specifications ..... 50

Identification label ..... 51

# WELLBOX® 2 WARRANTY CERTIFICATE

This warranty (the "Warranty") applies to all countries in which the Wellbox® 2 device (the "device") has been sold on the market by LPG Systems ("LPG") or an LPG® subsidiary, reseller and/or official distributor. LPG® warrants that the Device is, at the time of its original purchase, free from defects in material, design, and workmanship subject to the following provisions and conditions:

1. The Warranty is given to the purchaser. It in no way limits the purchaser's rights under the mandatory provisions of the applicable national law.

2. The warranty is valid as long as the purchaser provides original proof of purchase with the Device, including an invoice, receipt, sales slip, warranty certificate issued by LPG® along with the sales slip, or proof of payment by credit card if the purchase was made over the Internet. Proof of purchase must include the date of purchase, the specific name of the Device, and, if applicable, the serial number. LPG® reserves the right to refuse the free warranty if the above-mentioned documents cannot be presented or if the information they must contain is incomplete or illegible.

3. The warranty is only valid for private use of the device. It is excluded in the case of professional use or if the device has been rented on a commercial basis.

4. The legal warranty period is a maximum of two (2) years from the time the Wellbox® 2 is delivered to the purchaser. The warranty period may be extended by a commercial warranty of up to four (4) years if the purchaser chooses the extended warranty option when purchasing a new Wellbox® 2 from LPG®.

5. During the warranty period, LPG® will repair or replace the defective Device free of charge, at its sole discretion and while supplies last. LPG® will return the repaired device or another device in good working order to the purchaser. All parts or other materials that have been exchanged shall become the property of LPG®.

6. If, during the two years following purchase from LPG®, the device is repaired under the legal warranty of conformity, the buyer will benefit from an additional six (6) months of warranty. If the device is exchanged, the purchaser's legal warranty of conformity will be extended by two (2) years.

7. The Warranty applies to defects that appear when the device is used in accordance with the user manual. It does not cover maintenance or the repair or replacement of parts as a result of normal wear and tear or ageing of the device. Furthermore, the Warranty does not apply in particular when :

a. the defect was caused by the device being subjected to: use contrary to the user manual, rough handling, exposure to excessive humidity, extreme thermal or environmental conditions or rapid changes in these conditions, exposure to corrosion or oxidation, unauthorised modifications or connections, repairs carried out by unauthorised repairers or by the purchaser himself, repairs using unauthorised spare parts, misuse, unsuitable power supply, failure to follow assembly instructions, incorrect installation, accident, natural forces, contact with liquids, the action of chemicals and any other action reasonably beyond the control of LPG®. Use only on clean, dry skin. Do not use liquids and/or creams with Wellbox® 2.

b. LPG® was not notified by the buyer of the defect within thirty (30) days of the defect appearing during the warranty period;

c. The Product was not returned to LPG® within sixty (60) days following the appearance of the defect during the warranty period;

d. The device's serial number or any other identification number has been removed, erased, damaged, altered, or is no longer legible;

e. The defect results from the fact that the device has been used with an accessory not manufactured or supplied by LPG®.

8. The purchaser undertakes to return the device in its original packaging. LPG® assumes no liability in the event of damage to the Device due to transport in inappropriate packaging. LPG® reserves the right to request a flat-rate contribution from the buyer to cover product shipping costs.

9. The Warranty constitutes the purchaser's sole and exclusive remedy against LPG® and LPG®'s sole and exclusive liability to the purchaser for defects in the device that appear when the device is used in accordance with the user manual, as defined in Article 7.

The Warranty replaces any other warranty or liability, whether oral, written, legal (insofar as it is not mandatory), contractual, quasi-default or other. LPG® shall not be liable for any incidental, consequential or indirect damages or expenses. Any modification of the warranty requires the prior written consent of LPG®.

10. In accordance with the provisions of Article L111-4 of the French Consumer Code, the buyer is informed that spare parts essential to the use of the Wellbox® 2 device are available on the market for a period of three (3) years after the last new device is placed on the market.

LPG is registered as a producer under number FR002165\_05SRUY.

11. In the event that the purchaser must return the device to LPG® as part of the Warranty or for after-sales service, the purchaser agrees to erase all data contained on the device before sending the device to LPG® for repair.

Upon receipt of the Device at LPG®'s facilities, the Device will be reinitialized and will no longer contain the data initially entered by the purchaser when it is returned to LPG®.

You can activate your 2-year warranty on the [www.wellbox.com](http://www.wellbox.com) website or by scanning the QR code opposite using your phone's camera.

## PACKAGE CONTENTS

The Wellbox® 2 device comes with all the equipment you need for body and facial treatments.

- 3 ROLL care heads (1,2,3)
- 2 LIFT care heads (4,5)
- A connection hose for connecting the care heads
- A power supply unit
- Two disassembly tools for the LIFT 4 and 5 heads
- Replacement filters
- Replacement flaps for the LIFT heads
- A quick user guide

**IMPORTANT: Keep the original Wellbox® 2 packaging for the duration of the warranty. Any returns during this period must include the original packaging.**

## CONTRAINDICATIONS

Do not use the Wellbox® 2 device if you have any medical contraindications, or if you have not read this chapter carefully and reviewed the precautions in this guide. Consult your doctor if you have any doubts about the compatibility of your health with the use of the Wellbox® 2 device.

The Wellbox® 2 device can be used on all healthy external parts of the body, with the exception of the eyes, ears, glottis and genitals.

The Wellbox® 2 device should not be used on mucous membranes, nipples, the spine or on any suspicious moles/beauty spots without medical advice.

If you are pregnant, do not use the device on your stomach/abdomen or lower back. Some medical conditions may require the advice of your doctor.

Do not use the Wellbox® 2 device if you have any of the following:

- Epilepsy problems
- Blood clotting disorders
- Heart problems
- Pain or redness
- Tingling in the limbs
- Phlebitis, varicose veins or if you have undergone an operation on your veins in the previous two months.
- Fever, infection, inflammation, acute dermatitis, open wounds, recent scars or any serious and/or progressive illness.
- In case of any prolonged use of medications (anti-inflammatories, anticoagulants, antibiotics and antihistamines), consult your doctor before using this device.

Discontinue care if you experience any unexplained pain in your calves.

If you experience any side effects related to the Wellbox® 2 treatments during or after use, consult your doctor and contact Wellbox® Customer Service.

If you have any doubts about the compatibility of your health with the use of the Wellbox® 2 device, consult your doctor.

## PRECAUTIONS

- Do not use the Wellbox® 2 device on children under the age of 7. Do not allow a child to use the Wellbox® 2 device alone. Remember to reduce the intensity when treating a child.
- Close supervision is necessary when the device is used by, on or near children, invalids or disabled people.
- This device may be used by children aged 11 and over and by persons with reduced physical, sensory or mental capabilities or lack of experience and knowledge if they have been given supervision or instructions to use this device safely, and provided they understand the hazards involved.
- Do not allow children to play with the device.
- Cleaning and maintenance must not be carried out by unsupervised children.
- **Using the Wellbox® 2 device should never be painful or leave marks; any redness should not last more than a few seconds or leave traces.**
- Never use the Wellbox® 2 device without its protective filter.
- Use this device only as described in this guide. Use only the care heads and accessories supplied by LPG SYSTEMS.
- Do not use detergents, solvents, chemicals or abrasive sponges to clean the Wellbox® 2 device. Use wipes soaked with a bactericide and fungicide solution.
- Never use the Wellbox® 2 device on jewellery or piercings.
- Do not operate the device at altitudes above 2,000 metres.
- Under normal conditions of use, the device should be stopped for 30 minutes between two 30-minute sessions.
- Do not use the device in an oxygen-rich environment.
- The device should always be positioned so that the mains plug is easily accessible.
- **The device must always be used on clean, healthy skin (no creams, oils, etc. before use).**

**Warning: It is strictly forbidden to modify your device without the authorisation of LPG Systems**

Using the Wellbox® 2 device regularly will help you feel better; we recommend not interrupting its use for too long. Failing to follow the instructions and recommendations may prevent you from obtaining the desired results. Do not use the Wellbox® 2 device for longer than the recommended time. **Results may vary depending on your age, your lifestyle (sporting activity, diet, etc.) and your morphology.**

## IMPORTANT SAFETY INFORMATION

Before installation and use, you are strongly advised to read the instructions in this chapter.

### **DANGER - TO REDUCE THE RISK OF ELECTRIC SHOCK:**

- Check that the mains voltage is consistent with the range indicated for the device's power supply (100-240 VAC). The device must be connected to a wall socket that complies with your country's specific standards. Avoid connecting the Wellbox® 2 device to the mains via a power strip. Disconnect the Wellbox® 2 device from the power supply after each use or before cleaning.
- Do not use the Wellbox® 2 device while taking a bath or shower. Do not place the Wellbox® 2 device on a damp or unstable surface or near to a bathtub, shower or sink. Use the Wellbox® 2 device in a dry place. Avoid splashing any water on the Wellbox® 2 device. Do not touch any device that has fallen into water; disconnect the power supply immediately. Do not immerse the hose in water. Do not use the Wellbox® 2 device outdoors.
- Do not drop or insert objects into the openings of the device.
- Do not allow any solid object or liquid to be sucked into the Wellbox® 2 device.
- The device must always be used only with the power supply unit provided.
- The device must be powered only with the safety extra-low voltage corresponding to the marking on the device.

**WARNING - TO REDUCE THE RISK OF BURNS, FIRE, ELECTRIC SHOCK OR INJURY:**

- Do not use the device under a blanket or similar item. Excessive heat could occur and cause fire, electric shock or injury to persons.
- Never use the device if it has a damaged cable or plug, if it is not working properly, if it has been dropped or damaged, or if it has fallen into water. In these cases, contact Wellbox® Customer Service.
- Do not carry or move the device by its top cover, cables or hose, or use them as carrying handles. Do not drop the Wellbox® 2 device.
- Do not shake or knock the device.
- Do not use the Wellbox® 2 device without first connecting the hose to the heads. Keep cables and hoses away from sources of heat. Do not use the device when the ambient temperature is higher than 40°C.
- Do not use the device after using aerosol products (sprays) or when oxygen is being administered.
- Do not leave the device unattended while plugged in. To switch off the device, press the on/off button, then disconnect the AC/DC power supply unit from the mains socket.
- Do not open or disassemble the device.
- The Wellbox® 2 device must be powered only with 24V DC safety extra-low voltage.

**DESCRIPTION OF THE DEVICE**

1. Touch screen.
2. Removable hose for connecting the device to one of the care heads.
3. Connecting the hose to the device.
4. Storage for the four heads: ROLL heads 2 and 3, LIFT heads 4 and 5.
5. Lid release button.
6. Motorised ROLL 1 care head.
7. Hose storage compartment.
8. Protective filter access plug.
9. Protective filter accessible under plug (8).
10. Two removal tools for LIFT heads.
11. Carrying handle.
12. Connection point for power cable.

## DESCRIPTION OF THE TOUCH SCREEN

The Wellbox® 2 device is equipped with a touch screen with access to 5 menus:

**MY ACTIVITY:** find your current programmes, your care history, and track your progress.

**TARGETED CARE:** work on the area (body or face) of your choice depending on your objective (cellulite, slimming, radiance, firmness, etc.)

**PROGRAMMES:** follow a programme for several weeks (face, body or body & face)

**SETTINGS:** define your settings and get information on your device

**PROFILE:** modify your profile or create a new one (maximum two profiles)

For step-by-step guidance, see the full user guide available on our website.

Watch video advice from the Wellbox® 2 expert by scanning the QR code opposite using the camera of your phone

## DESCRIPTION OF CARE HEADS

The Wellbox® 2 is fully equipped with 5 care heads (3 ROLL heads and 2 LIFT heads) so that you always have the solution best suited to the desired objective and the area to be treated.

3 Roll heads ① ② ③

Each Roll head has removable sealing blades.


Head Roll ①

1. Trigger for reversing the direction of the motorised rollers indicated by an arrow.
2. Motorised rollers.
3. Treatment chamber.
4. Removable sealing stripes (behind protective parts).

This motorised Roll 1 head can treat all areas of the body.

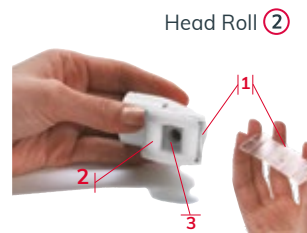Head Roll ②

1. Removable sealing strips.
2. Removable rollers.
3. Treatment chamber.

This small Roll 2 head can be used for fatty areas of the face (double chin).

Head Roll ③

1. Removable sealing strips.
2. Removable rollers.
3. Locking button.  
When the button is pressed ←, the rollers are locked (continuous mode), when released →, the rollers are free to move towards or away from each other (rhythmicity).
4. Treatment chamber.

The Roll 3 head can be used on the arms, knees or calves depending on body shape.

2 Lift heads ④ ⑤

Lift head ④

Lift head ⑤

- The LIFT 4 head is used for large areas: forehead, cheeks, oval, neck, décolleté.
- The LIFT 5 head is used for narrower areas: eye contour, frown lines (between the eyebrows), mouth contour, hands.

## FIRST INSTALLATION

- 1 Open the device by pressing the push-button on the front.
- 2 Store the 5 care heads in their dedicated slots, recognisable by their markings.
- 3 Insert the hose with the raised arrow pointing upwards as shown in the photo.
- 4 Choose one of the care heads (LIFT or ROLL) and connect it by aligning the red arrows.
- 5 The Wellbox® 2 device is supplied with a power supply unit. Input 100-240 VAC, 50/60Hz - Output 24 VDC. Connect the power cable to the back of the Wellbox® 2 device. Then plug the power supply unit into a wall socket.
- 6 Once you have plugged in your device, the screen will turn on automatically. Now you can enjoy your device!

For optimum performance, we recommend that you leave your device plugged in for at least 2 hours:

- After the first use
- After any long absence of use (6 months)

**Warning: The power supply unit is an integral part of your device. It is forbidden and dangerous to dismantle this box or to use a power supply other than that supplied with the device.**

## START-UP

- 1 Once the power supply has been connected, the device's touch screen will light up.
- 2 You will need to fill in a number of details when you first start up to create your profile: language, date, time, first name, gender, screen colour. You are the sole owner of the rights to the data processed in connection with the use of the Wellbox® 2 device (hereinafter referred to as the "Data"). Access to the Data is reserved for you. Neither LPG® nor its distributors/resellers have any access to the Data collected, processed, and/or stored on the Wellbox the Wellbox® 2 device.

You are solely responsible for the use of the Wellbox® 2 and, in particular the quality, lawfulness and relevance of the Data it processes.

You are responsible for backing up your device. LPG® cannot be held responsible for any loss and/or alteration of your device. If your Wellbox® 2 is returned to LPG® for any reason, it will be automatically reset (resulting in the deletion of all of the device it contains) upon receipt at LPG®'s premises.

- 3 The home page features 4 menus (MY ACTIVITY, TARGET CARE, PROGRAMS, SETTINGS). Depending on your objective, you can choose between a targeted care or a programme over several weeks.
- 4 Once you have chosen your care or programme, the screen will display the care head(s) to be used. Plug in the corresponding care head and press PLAY to start the care.
- 5 You can increase or decrease the intensity of the care according to how you feel and the quality of your skin by pressing + or -. Note that the care must never hurt.
- 6 The care steps will follow automatically. An audible beep informs you to the next step.
- 7 You can stop the care at any time and resume it whenever you like by pressing PAUSE.

- 8 When you have finished using the device, press the OFF button at the top left of the screen to switch it off and unplug the power cable. If the power cable is not unplugged, the device will go into standby mode automatically. To return to the home page, press the screen.

Watch video advice from the Wellbox® 2 expert by scanning the QR code opposite using the camera of your phone

## STORAGE

To close the Wellbox® 2 device after use:

- 1 Connect the ROLL 1 head to the hose.
- 2 Wind the hose into the purpose-built groove (2 turns).
- 3 Then place the ROLL 1 head in the space reserved for it next to the touch screen.

In case of any difficulty, we recommend removing the ROLL 1 head from the hose or removing the hose completely from the Wellbox® 2 device.

1

2

3

Watch video advice from the Wellbox® 2 expert

# MAINTENANCE

For reasons of hygiene, but also to keep it in good working order and ensure it lasts as long as possible, you need to clean your device regularly.

**IMPORTANT: The Wellbox® 2 must be disconnected from the mains before any cleaning operation.**

## 1. CHANGING THE PROTECTIVE FILTER

The Wellbox® 2 has a protective filter to ensure that it works properly. When this filter no longer allows you to use your device correctly, an indication appears on the screen to replace it.

To change the filter :

- 1 Unplug the power supply from the wall socket.
- 2 Unscrew the plug 1/4 turn to the left as shown in the photo then pull it towards you.
- 3 Remove the old filter by loosening it from the plug and discard it.
- 4 Refit the new filter (only one at a time) by folding back the 2 lugs then replace the plug.

Watch the video advice on how to maintain your Wellbox® 2 by scanning the QR code opposite using your phone's camera.

## 2. ROLL HEAD MAINTENANCE

For reasons of hygiene, the care heads must be cleaned before and after each use, with the heads disconnected, using wipes soaked with a bactericide and fungicide solution.

To clean the Roll head ①

- 1 Holding the head with one hand, unclip the protective part from each roller, then wipe them clean.
- 2 Then remove the 2 sealing strips by carefully gently pulling them out of their notches, then clean them carefully with a cloth.
- 3 Clean the treatment chamber and rollers with a wipe.
- 4 Clip the sealing strips onto the protective parts with the forked end on the inside. Insert the assembly into the head. Apply pressure on the notches to snap the assembly back into place.

To clean the Roll heads ② ③

- 1 Hold the head with the rollers facing you. Remove the rollers by pressing outwards.
- 2 Remove the sealing stripe from the treatment chamber.
- 3 Clean the treatment chamber, the sealing strip and the rollers using a wipe.
- 4 Replace the sealing strip in its original position.
- 5 Replace the rollers, pressing inwards until they click into place.


### 3. LIFT HEAD MAINTENANCE

For reasons of hygiene, the care heads must be cleaned before and after each use, with the heads disconnected, using wipes soaked with a bactericide and fungicide solution.

To clean the Lift heads ④ ⑤

- 1 Take the head Lift ④ or ⑤ and insert the tool [1] ou [2], specific to each head, push the bevelled side into the hole in the removable flap until it clicks.
- 2 Take out the tool; the flap comes with it.
- 3 Remove the flap from the tool and clean it thoroughly with a wipe.
- 4 Take the tool on the smooth side and replace the flap on the tool.
- 5 Insert it into the care head until it stops at the bottom of the head.

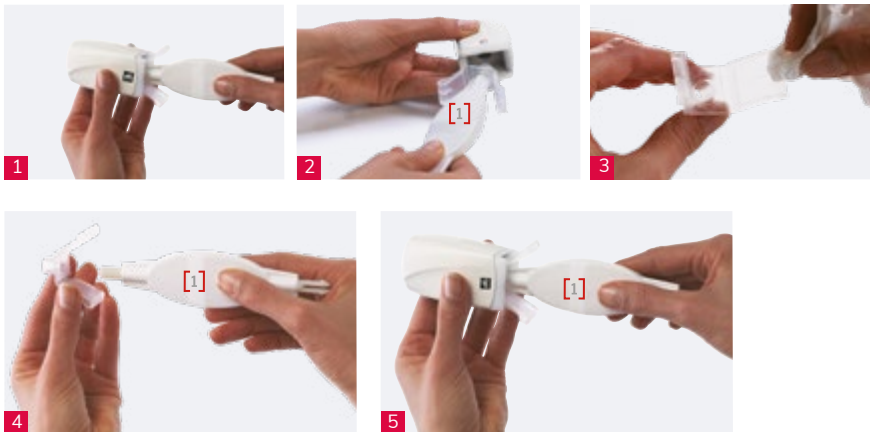Watch the video advice on how to maintain your Wellbox® 2 by scanning the QR code opposite using your phone's camera.

### A PROBLEM? WHAT TO DO?

If your device is not working properly, we recommend that you carry out the following checks before contacting Wellbox® 2 Customer Service:

- Have the precautions for use been observed?
- Is the device connected to the power supply?
- Is the power supply unit suitable for use with the device?
- Does the power cable you have selected comply with the standards in your country?
- Has the Wellbox® 2 been unused for five minutes?
- Does the screen light up?
- Is the filter clean and in place?
- Is the hose not blocked? Is it connected correctly?
- Is the care head correctly connected to the hose?
- Are the parts of the care heads (sealing stripes, rollers, flaps) correctly fitted, clean and mobile?
- Does the control button on the motorised Roll head 1 operate effortlessly?
- For safety reasons, do not dismantle either the power supply or the device. Otherwise, the warranty will not apply.

Once these checks have been carried out, if the malfunction persists, contact Wellbox® 2 Customer Service or an authorised distributor, quoting the device serial number shown on the barcode label on the underside of the device (for example: WB2EDL1244158) and its operating time in hours.

- [www.wellbox.com](http://www.wellbox.com)
- If you have any doubts about the operation or use of your device, please contact Wellbox® Customer Service: +33(0)4 75 78 69 00 (toll-free in Europe only)

**Warning: Please check that you have followed the installation and operating instructions beforehand and use instructions.**

## TECHNICAL SPECIFICATIONS

Dimensions L x W x H: 350 x 300 x 250 mm

Net weight: 7.5 kg

Power consumption: 70 Watts

Class 3 device, maximum duration of use per session: 30 minutes

Time required for the Wellbox® 2 device to warm up if stored at -25°C: 9 hours

Time required for the Wellbox® 2 device to cool down if stored at +70°C: 18 hours

Maximum surface temperature of the Roll 1 head: 42.5°C

Maximum surface temperature of the rollers on the Roll 1 head: 44.5°C

Protection rating \_\_\_\_\_ IP21

Electrical characteristics \_\_\_\_\_ 24V; 3.75A

Maximum regulated vacuum \_\_\_\_\_ 350mbar

Maximum altitude \_\_\_\_\_ 2.000 m

### TECHNICAL CHARACTERISTICS OF THE POWER SUPPLY UNIT

Manufacturer \_\_\_\_\_ MeanWell

Ref. \_\_\_\_\_ GSM90B24-P1M

Electrical protection class \_\_\_\_\_ II

Protection rating \_\_\_\_\_ IP22

Electrical specifications \_\_\_\_\_ 100-240VAC 50/60Hz 1.5A

### ENVIRONMENTAL CONDITIONS OF USE

Temperature \_\_\_\_\_ +5°C to +30°C

Relative humidity \_\_\_\_\_ 15% to 90%  
(without condensation)

Atmospheric pressure \_\_\_\_\_ 700 hPa to 1.060 hPa  
(for use in a normally ventilated room).

### ENVIRONMENTAL CHARACTERISTICS OF TRANSPORT AND STORAGE

Temperature \_\_\_\_\_ -20°C to +70°C

Relative humidity \_\_\_\_\_ 10% to 90%  
(without condensation)

Atmospheric pressure \_\_\_\_\_ 500 hPa to 1.060 hPa

### IDENTIFICATION LABEL

"Type" identification label subject to change to take account of upgrades to the device and to the Wellbox® 2 device type. The label on the bottom of your Wellbox® 2 device is the approved version.

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**WARNING**  
The manufacturer reserves the right to change the product specifications of the device.

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are trademarks of LPG Systems.
**SIÈGE SOCIAL LPG SYSTEMS SAS
LPG SYSTEMS SAS HEAD OFFICE**