wellbox S

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QUICKSTART

AVVIO RAPIDO

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SCHNELLSTART

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Ref. SG0906-ED-IndC3-062025

[www.wellbox.com](http://www.wellbox.com)

LPG

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Thank you for acquiring Wellbox®[S], a device for at-home use, whose effects have been scientifically demonstrated. Your Wellbox®[S] device will give you full satisfaction if you follow the protocols included in the user guide.

## USER GUIDE

- > To reduce the environmental impact while still maintaining the level of safety and ease of use, the complete user guide is available in electronic format at [www.wellbox.com](http://www.wellbox.com).
- > This guide must be read before you use your device. If you have the slightest doubt about the operation or use of your device, or if you wish to receive these instructions on paper, do not hesitate to contact


**Wellbox® Customer Service hotline +33(0)4 75 78 69 00**  
(toll-free in Europe only)

## PRODUCT WARRANTY

Also, take the time to activate your warranty directly on line.

[www.wellbox.com](http://www.wellbox.com)
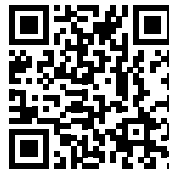
### Identification label

"Type" identification label subject to change to reflect changes to the device and the Wellbox's type. The label at the bottom of your Wellbox is the approved version.

### Description of symbols

	Alternate current		Fragile, handle with care
	Continuous current		Not reverse
	Warning; refer to the additional documentation		Store the device somewhere protected from inclement weather
	Standby, moving from standby mode to full power		Temperature limits
	Notice		Relative humidity limits
	Dangerous voltage		Recyclable packaging, it is the subject of a sorting order
	This icon indicates that the unit was sold after August 13, 2005.		Recyclable packaging
<p>In conformity with the 2002/96/CE electronic waste directive, your device includes components that cannot be thrown away with standard household waste but must be disposed of by means of recycling. At the end of the service life of your device, for proper recycling, return the device to your supplier or to a designated collection point. By doing so, you help the environment by contributing to the conservation of natural resources and the protection of human health.</p>			Consult the accompanying documents
			Indoor use only
			Power supply connector polarity
			Double isolated
			Class 3 device

The manufacturer reserves the right to change the technical specifications of the device.

*Read this Installation and User Guide carefully and keep it.*

## CONTRAINDICATIONS

Wellbox®[S] can be used on all healthy external parts of the body, with the exception of the eyes, ears, glottis, and genitals. The device should not be used on mucous membranes, nipples, the spine or on any suspicious moles without medical advice.

### Do not use Wellbox®[S] in case of:

- Epilepsy problems
- Blood clotting disorders
- Heart problems
- Pain or red blotches
- Tingling in the limbs
- In case of pregnancy, do not use the device on the areas of the abdomen and lower back. Some conditions may require the advice of your doctor.
- Phlebitis (blood clots), varicose veins, or if you have undergone an operation on the veins within the last 2 months
- Fever, infection, inflammation, acute dermatitis, open wounds, recent scars, or any serious and/or progressive illness
- If taking long-term medications (anti-inflammatory drugs, anticoagulants, antibiotics or antihistamines)
- Wellbox®[S] cannot be applied on areas that are swollen, or affected by rosacea

Stop the treatment in the event of unexplained calve pain. If any side effects linked to Wellbox®[S] treatments appear during or after use, consult your doctor and contact LPG® Customer Service Department.

If you have the slightest doubt concerning the compatibility of your state of health with the use of Wellbox®[S], seek your doctor's advice.

## FIRST USE

### > TREATMENT HEADS:

WATCH THE VIDEO!

Open the device by pressing on the button located in front. Store the four treatment heads in the locations specifically designed for that purpose. Each head's dedicated storage area is indicated with a pictogram.

### > INSTALLATION:

Insert the hose with the arrow  (in relief) upwards, as indicated in the photo.

Choose one of the heads ( LIFT or ROLL) and then connect it to the hose by lining up the arrows on the head and hose.

### > POWER SUPPLY:

Insert the male jack into the female socket located at the rear of the Wellbox®[S device]. Connect the power supply to a wall outlet by first choosing the right one for your country.

### > STARTING UP:

After connecting the power, press the key  to turn on the device. The control panel lights up.

The Wellbox®[S] Routine, selected by the treatment menu will inform you, all the power settings to adjust.

- The suction power (Power) must be adjusted depending on how you feel and on your skin quality.
- Suction sequentiality (Select) is the number of aspirations per second.

- continue suction
- 4 micro-aspirations per seconde
- 8 micro-aspirations per seconde
- 16 micro-aspirations per seconde

These two pictograms ROLL head and LIFT head are provided for information purposes only. They will guide you through the routine and indicate which treatment head to use according to the suction sequentiality you have selected.

• Treatments:

Select the treatment heads (LIFT or ROLL) to use, and then adjust the power settings (suction power, suction sequentiality, The duration of the treatment) according to the treatment menu.

FACE TREATMENT

by the white arrow.

For better stimulation, move your LIFT head with small hops, stopping 4 seconds on each point, while respecting the direction indicated


The red lines indicate how to position the flaps when moving the treatment head. They must always be perpendicular to the wrinkle.

BODY TREATMENT

The white arrows indicate the direction that the treatment head must follow. When using the ROLL(3) head, use the trigger to reverse the direction of the rollers.


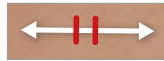
The red lines indicate how to position the rollers when moving the treatment head.

Press to start the treatment. You can interrupt the treatment at any time and resume when you want by pressing .



> BEAUTY ROUTINES

To meet all your beauty objectives, Wellbox® [S] offers a multitude of body and face routines. Among them:

BODY ROUTINE : FLAT STOMACH

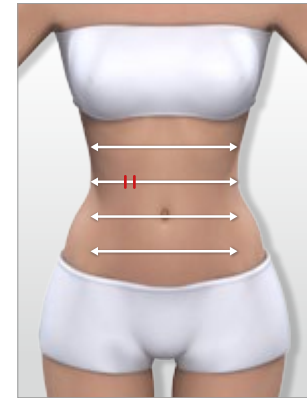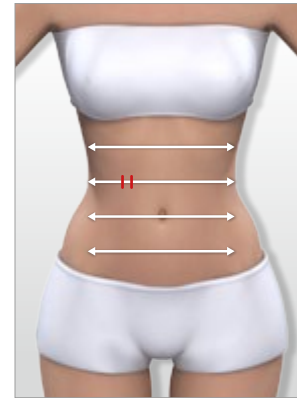FACE ROUTINE : FACE FIRING

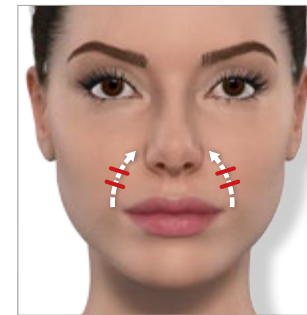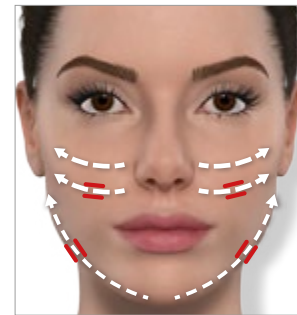Find all the beauty routines here:

> STORAGE:

To close your Wellbox®[S] device after use:

- Connect the ROLL(3) head to the hose.
- Roll the hose into the tube designed for this purpose (two times around).
- Place the ROLL(3) head in the location next to the control panel.

## IMPORTANT SAFETY INSTRUCTIONS

- > Check that the network voltage complies with the range indicated on the power supply of the device (100-240 VAC). It should be connected to a wall socket that complies with the electrical standards in force. Avoid connecting Wellbox®[S] to a multi-socket power bar. After every use and before cleaning, disconnect the Wellbox®[S] power system.
- > Do not use Wellbox®[S] while taking a bath or shower. Do not place Wellbox®[S] on a wet or unstable surface, or near a bathtub, shower or sink. Use Wellbox®[S] in a dry area. Avoid splashing any water on your device. Do not touch a device that may have fallen into water. Immediately disconnect from the power supply. Do not put the hose in water. Do not use Wellbox®[S] outdoors.
- > Do not drop the device, or insert objects in any openings. Do not suck up any solid bodies or liquids with Wellbox®[S].

### **WARNING: In order to reduce any risk of burns, fire, electric shock or injuries:**

- > Do not use the device under any kind of blanket or cover. Excessive heat may be produced which may provoke fire or an electric shock, or injure someone.
- > Never use the device if a cable or connection is damaged, if it doesn't work properly, if it has been dropped or damaged, or if it has fallen into water. In these cases, please contact LPG® Customer Service.
- > Do not pick up or move the device by its cover, cables or hose, or use them as carrying handles. Do not drop Wellbox®[S]. Do not shake or hit the device.
- > Do not use Wellbox®[S] without first connecting the hose to the heads. Keep the cables and the hose away from any source of heat. Do not use the device when the temperature is above 40°C.
- > Do not use the device when aerosol sprays have been used or when oxygen has been administered.
- > The device must not be left unsupervised when connected. To turn the device off, press the on/off push-button, then disconnect the AC/DC power supply from the wall socket.
- > Do not open or dismantle the device.

## SAFETY PRECAUTIONS

- > Do not use Wellbox®[S] on a child under the age of 7, and do not let a child use Wellbox®[S] alone. Remember to reduce the intensity when treating a child.
- > Close supervision is necessary when the device is used by, on or near children, or disabled or handicapped persons.
- > Using Wellbox®[S] should never be painful or leave marks; any redness should disappear after a few seconds and there should be no marks.
- > Never use Wellbox®[S] without its protection filter.
- > Use the device only as described in this guide. Use only the treatment heads and accessories supplied by LPG SYSTEMS.
- > This device is exclusively for at-home use.
- > Do not use detergents, solvents, chemical products or abrasive sponges to clean Wellbox®[S]. Use wipes soaked in a bactericidal and fungicidal solution.
- > Never use Wellbox®[S] on jewellery or piercings.
- > Do not operate the device at an altitude above 2,500 m.
- > In normal use, the device must be stopped for 30 minutes between each 30-minute session. Using Wellbox®[S] regularly makes you feel better; it is recommended that its use not be interrupted for too long a period. If the instructions and recommendations for use are not followed, you risk not achieving your desired results. Do not use Wellbox®[S] for longer than the recommended time.

Results can vary according to age, lifestyle (sports activity, diet, etc.) and morphology.

**Caution: It is strictly forbidden to modify your device without the permission of LPG SYSTEMS.**

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**Siège social LPG SYSTEMS SAS**